

Thinking

NAUCO  
**inserto**<sup>®</sup>  
SEAT




CERTIFIED COMPANY WITH QUALITY  
MANAGEMENT SYSTEM ISO 9001 / EN ISO 13465  
N° 279412 / N° 279411

2 A Teschio e denti, 1489.

RL 19058 v (B. 41 r; K/P 42 v).

Il vacuo della cassa dell'occhio, e 'i vacuo dell'osso sostenitore della guancia, e quello del naso e della bocca sono d'eguale profondità e terminano sotto il senso comune [sede centrale delle sensazioni] per linea perpendicolare. E ciascuna d'esse vacuità ha tanto di profondità quant'è la terza parte del volto dell'omo, cioè dal mento a capegli.

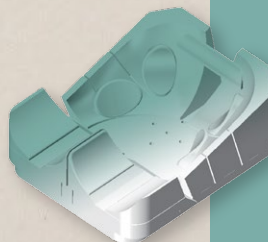

SEAT

A NEW SOLUTION  
INDIVIDUALLY CUSTOMISABLE

# A NEW SOLUTION

## INDIVIDUALLY CUSTOMISABLE

Construction kit of the individually customisable seat  
Perfect reproduction of anatomical shapes



### Why Nauco ?

- SUPPORTIVE
- LIGHTWEIGHT
- MULTIFUNCTIONAL



MODULAR AND CONFIGURABLE STRUCTURAL KIT



Free shapes PADDING KIT



STABILITY AND POSTURAL CONTROL, EXCEPTIONAL FEMORAL CONTROL

### CLINICAL INFORMATION / DESTINATION OF USE

The pelvis positioning system Inserto Nauco is recommended preferably, but not exhaustively, for users with mobility limits and poor postural control, in particular:

- those who present leg length discrepancies, wind-swept, pelvic obliquity, pelvis rotation, anterior and posterior pelvis tilt and risk of hip dislocation;
- those who need lateral contouring control of the legs and lumbar/sacral support in order to improve the self control of the trunk;
- for even very small users who have to be led to a sitting position;
- for people using the orthosis for the rachis;
- for users presenting a sacral/lumbar protrusion.

**PTS Pelvic Total Support**

The pelvis positioning system has to be used with any supporting / wheelchair base properly fixed through suitable gripping systems.

### TECHNICAL FEATURES

Each positioning system Inserto Nauco includes:

- a structural kit with different density combination composed by a flat base of construction and 22/24 polyethylene closed celled inserts, customisable and modifiable, useful to individually support or correct the posture of the pelvis during the life time of use;
- padding kit, washable at 60° with light detergent, which is composed by different layers:
  - a memory foam layer within the cover;
  - a flexible expanded polyurethane foam layer;
- Fireproof and antibacterial cover with an insider layer which makes it incontinent. It can be washed in tepid water with a light detergent. Do not tumble dry. The fabric doesn't require the use of disinfectants.

#### RICAMBI

- Flat base of construction
- Cover
- Padding kit

### WARNINGS


The Pelvis Positioning System Inserto Seat Nauco, in its design integrity and numerical and dimensional completeness of the supplied components (structural kit, padding kit, cover) can be easily adapted to sizes/morphology/deformities of the user. This kind of operations make the commissioning referable to a serial device.

Alternatively

The commissioning of the Positioning System Inserto Seat Nauco, deprived or modified, even partially, of its design integrity of numerical and dimensional completeness of the supplied components, built as per written prescription of a professional operator in function of the anatomy / morphology / deformities of the user, through the measurement detection and direct trials, can be referable to a custom device.

For partially sighted users this document can be visualised in PDF format (with the possibility to enlarge the text) on the web site [www.versainser.to.com](http://www.versainser.to.com)



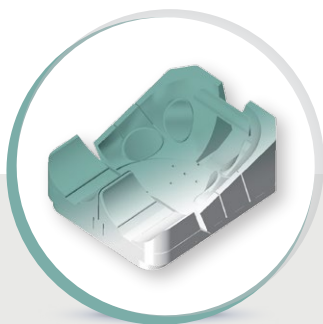

# Inserto Nauco

22/24\* MODULAR AND CONFIGURABLE INSERTS

COVER WITH AN INSIDER LAYER WHICH MAKES IT INCONTINENT


° Each postural kit is ordered based upon the size of mobility / wheelchair device.

INSERTO NAUCO SIZES TABLE	PAEDIATRIC				SIZES (cm) ADULT
MODEL	2225	2530	3034	3638	4042
Effective width (cm)	22	25	30	36	40
Max width achieved with adaptation (cm)	24	27	32	38	42
Min. width achieved with removal of material (cm)	20	23	28	34	38
Effective Length (cm)	25	30	34	38	42
Max. length achieved with adaptation (cm)	27,5	32,5	36,5	40,5	44,5
Min. length achievable	22	27	31	35	39
Min. and Max. range of anterior height of flat base of construction (14 cm)**	4/5 ** height referred to models 2225 and 4042				
Min. and Max. range of posterior height of flat base of construction +buttock inserts (9 cm)**	11,5/18,5 ** height referred to models 2225 and 4042				

\* possible and further reductions due to removal of material

◦ Min/Max positioning system weight: 0,4 kg/1,3 kg  
◦ Max. load referred to model 4042: 85 kg

conversion table	inches	6	7	8	9,5	10,5	11	12	12,5	14	16
	cm	15	18	21	24	27	28	30	32	36	40

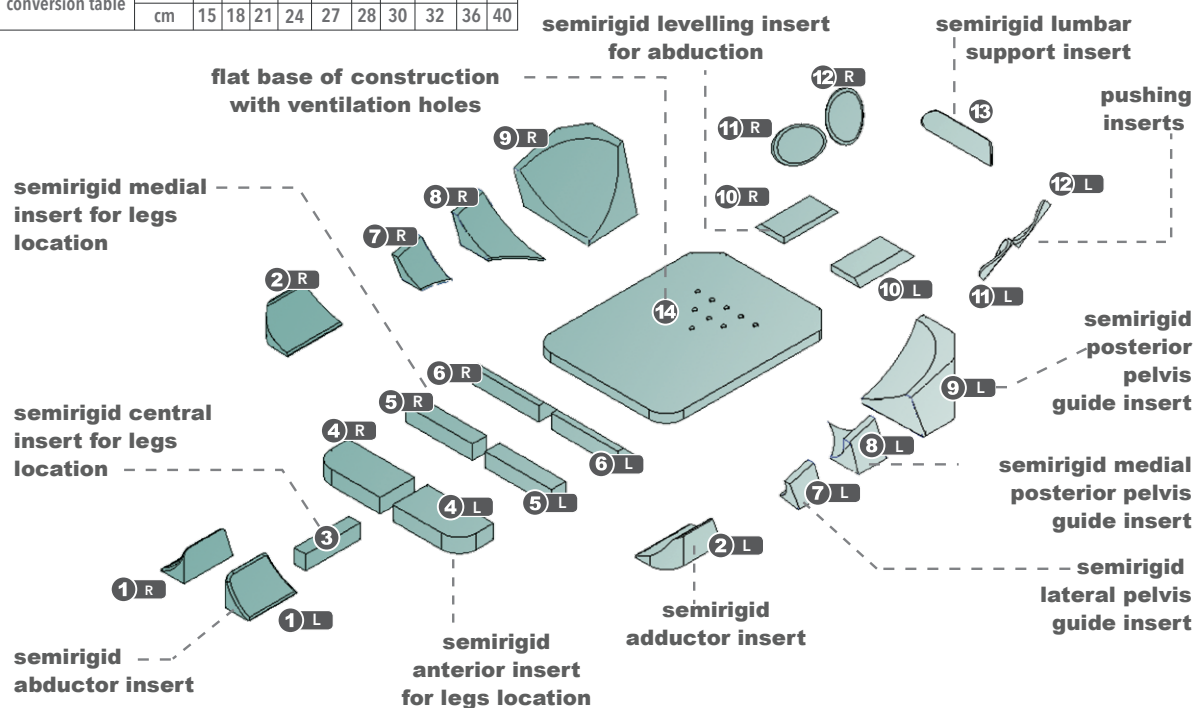




NAUCO




## Customisation of Nauco



HIP DISLOCATION PREVENTION

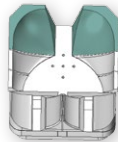

LOCATION FOR THE PROTRUSION

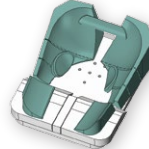


ADDUCTION




ABDUCTION




WIND SWEPT

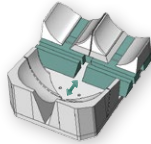

ANTI-THRUST OF THE PELVIS

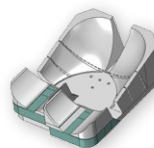


LOWER LIMBS ASIMMETRY

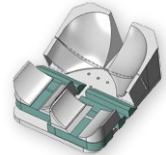


ASIMMETRY FOR HIP DISLOCATION OUTCOMES




TRUNK EXTENSION STIMULATION

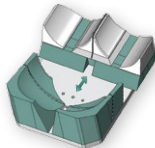


ANTERIOR-POSTERIOR PELVIS TILT

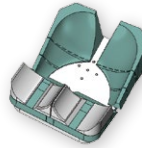


PELVIC ROTATION

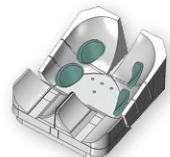

INCREMENT OF PELVIS CONTAINMENT

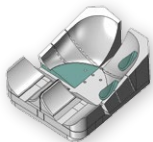

PELVIC OBLIQUITY

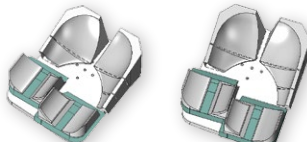

ASIMMETRY FOR PELVIS ROTATION

INDIVIDUALLY CUSTOMISABLE

Construction Kit of the individually customisable seat



# THE INNOVATIVE SOLUTION

A range developed to meet the individual postural needs (as per prescription) and numerous pathologies throughout the whole life cycle of use.

In order to maintain its features it is necessary to periodically check the conservation and wear of all the parts composing the device as well as all hook and loops stripes.

The positioning system seen the sizes and weight results easily transportable.

[www.versainsero.com](http://www.versainsero.com)

## 1 MAN/WOMAN ADAPTABILITY

Finally a solution that makes the difference with a fully personalised fit.

## 2 ADAPTABILITY TO THE INDIVIDUAL MORPHOLOGY AND ANATOMY

Unlimited possibilities to realize the customised postural solution.


## 3 PTS PELVIC TOTAL SUPPORT

A perfect fit to anatomy of the pelvis and related body segments, encouraging an active response of the musculo-skeletal system.

The distribution of the loads provided by the PTS reduces the mechanical stress on the deep tissues and significantly aids the prevention of decubitus ulcers.

## 4 FREE SHAPES PADDING

The innovative padding ensures maximum comfort and reliable protection.

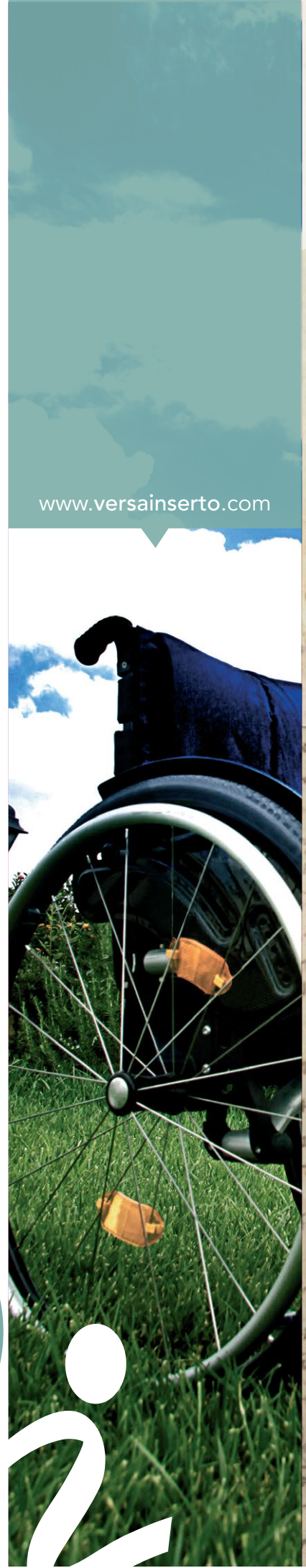



MANUFACTURED BY



Pro Medicare srl

Via Montagna Z.I. - Lotto 41  
72023 Mesagne (BR) - ITALY

Tel. +(39) 0831 777840

Fax. +(39) 0831 730739

Email: sales@promedicare.it


www.versainsero.com

www.promedicare.eu


CERTIFIED COMPANY WITH QUALITY  
MANAGEMENT SYSTEM ISO 9001 / EN ISO 13485  
N° 279412 / N° 279411



ISO 9001  
ISO 13485  
BUREAU VERITAS  
Certification

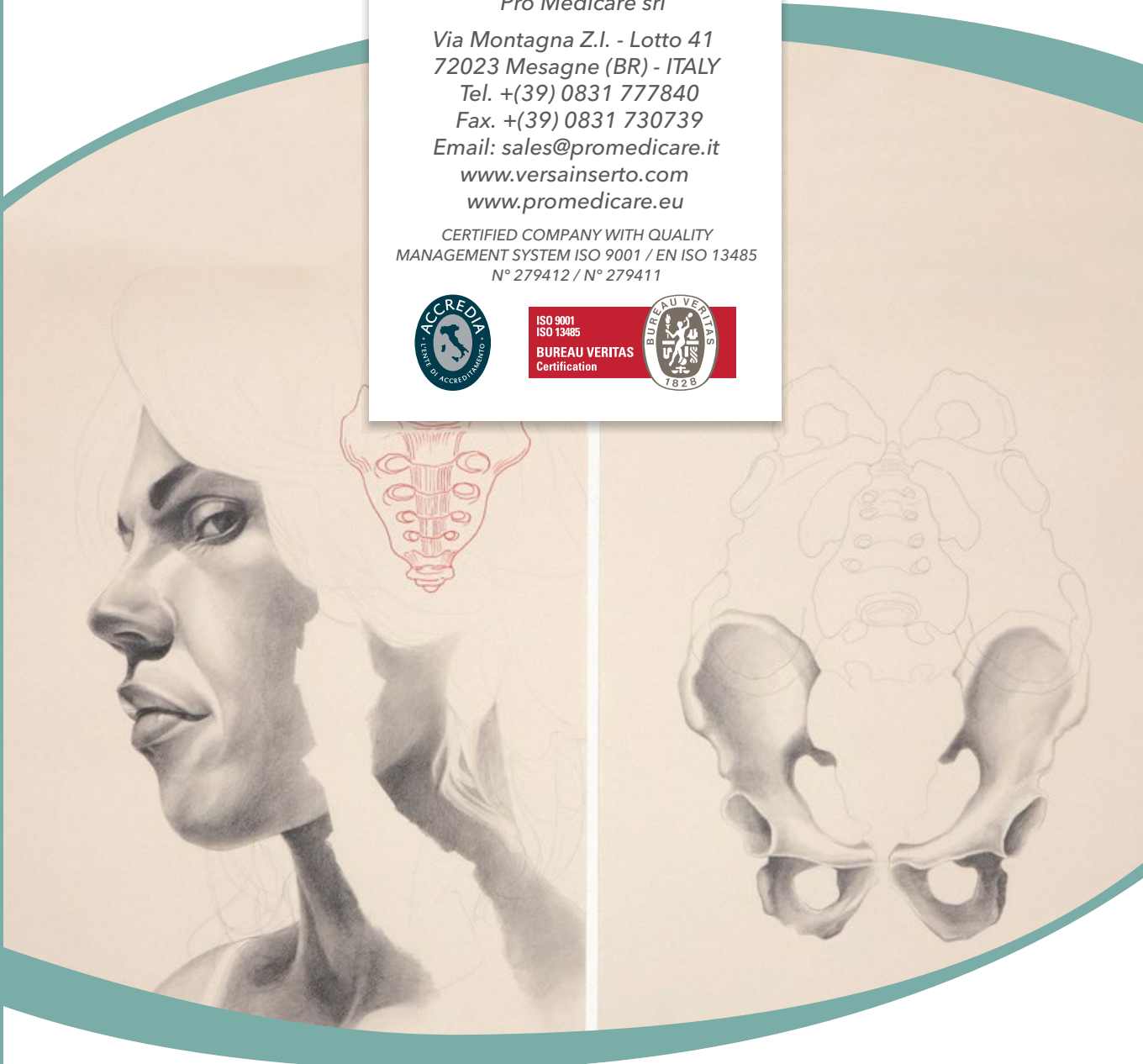

DISTRIBUTED BY

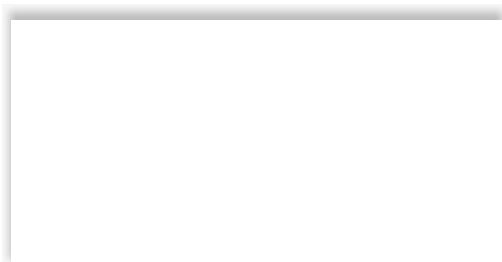

CE 2 YEAR WARRANTY

www.versainsero.com